



Professional

The Professional Ophthalmoscope has been designed to ensure ultimate control. Lenses, graticules and filters can all be introduced at the touch of a fingertip without the need to pull away from the patient.

- **Positive action maywheel and comprehensive lens range**

The positive action, easy access, maywheel not only allows you to dial up the required lens quickly, but also gives you total control of where you are in the lens range making each exam quicker and easier.

The extensive range of +29D to -30D in single dioptre steps ensures easy examination of all ocular structures and covers a wide range of patient/user prescriptions. A quick step swing of +/- 20D helps speed up diagnosis.

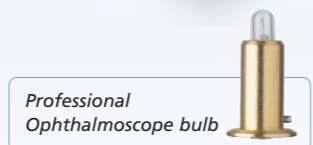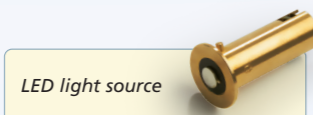

Key features

- **Choice of illumination**
The Professional now has a choice of Xenon or LED illumination for a brighter and whiter long lasting light. Giving you 3 x longer battery life.
- **Swing over filter**
The red free filter for enhanced vessel examination can be used in conjunction with any of the 6 diagnostic beams.
- **Brow rest**
The integral soft brow rest allows for easy positioning and protects your lenses if you are a spectacle wearer.



Part Numbers

Professional Ophthalmoscopes

1128-P-1014	3.6V Professional Ophth in case
1128-P-1015	3.6V LED Professional Ophth in case
1128-P-1013	2.8V Professional Ophth in case
1128-P-1008	2.8V LED Professional Ophth in case

Accessories




1011-P-7114	Std/Pract/Prof Xenon Bulb 3.6V Pkt2
1011-P-7106	Std/Pract/Prof Xenon Bulb 2.8V Pkt2
1011-P-7210	Std/Pract/Prof LED Bulb 3.6/2.8V Pkt1



Professional Retinoscope

Optimised for outstanding performance with either a Streak or Spot halogen beam with a simple bulb change. Plus the LED Professional Streak Retinoscope is now available!

- **Neutralisation check (Patented)**
The only Retinoscope allowing rapid confirmation of neutralisation. You can change the beam from divergent to convergent enabling you to cross check the accuracy of the neutral point, guaranteeing the assessment is correct every time.
- **Precision magnetic control (Patented)**
Improved 360° access to a single control for both rotation and convergence adjustments. The precision movements are achieved by use of rare earth magnet technology.

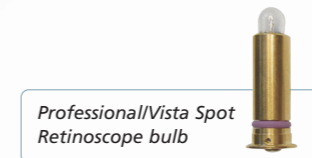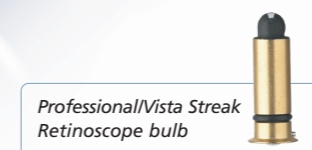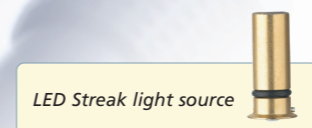

Key features

- **Choice of illumination**
Choose between 600 lux of brilliant halogen illumination or the LED Retinoscope. Whichever you choose, the Professional guarantees you enough light to refract even the most difficult eyes.
- **LED Professional Streak Retinoscope**
Gives you better contrast, image quality and longer battery life.
- **Two position aperture**
A simple finger movement allows a choice of 4mm or 1.7mm aperture to optimise brightness or resolution of the retinal reflex.

Part Numbers

Professional Retinoscopes

1302-P-1014	3.6V Professional Streak LED Ret and case
1302-P-1011	3.6V Professional Streak Retinoscope in case
1302-P-1009	3.6V Professional Spot Retinoscope in case
1302-P-1013	2.8V Professional Streak LED Ret and case
1302-P-1010	2.8V Professional Streak Retinoscope in case
1302-P-1008	2.8v Professional Spot Retinoscope in case

Accessory

1305-P-7000	Streak Retinoscope LED blister pack
-------------	-------------------------------------

Professional Combi Retinoscope Bulbs

1013-P-7009	Vista/Professional Streak bulb 3.6V Pkt2
1013-P-7011	Professional Spot bulb 3.6V Pkt2
1013-P-7010	Professional Spot bulb 2.8V Pkt2

Vista Retinoscope Bulbs

1013-P-7009	Vista/Professional Streak bulb 3.6V Pkt2
1013-P-7007	Vista Spot halogen bulb 3.6V Pkt2
1013-P-7006	Vista Spot halogen bulb 2.8V Pkt2